



What can I do if I have diabetes?

A. Use medicines wisely

- Sometimes people with diabetes need to take pills or insulin shots. Follow the directions.
- Ask your doctor, nurse, or pharmacist what your medicines do. Also ask when to take them and if they have any side effects.

B. Watch what you eat and get exercise

- There is no one diet for people with diabetes. Work with your health care team to come up with a plan for you.
- Be active at least 30 minutes a day, most days of the week. Exercise helps your body use insulin better.

Check your blood sugar and know your ABC's

- Help prevent heart disease and stroke by keeping your blood sugar, blood pressure, and cholesterol under control.
- Check your blood sugar with a test you can use at home.
- Ask your doctor for an A-1-C ("A-one-see") blood test. It checks blood sugar levels over 2 to 3 months.
- Talk to your health care team about your ABC's:
A-1-C
Blood pressure
Cholesterol

Reference: Food and Drug Administration (FDA) Diabetes Information

www.bionime.com

BIONIME

BIONIME CORPORATION

No.694, Renhua Rd., Dali Dist.,
Taichung City 412, Taiwan (R.O.C.)
Tel: +886 4 24951268
Fax: +886 4 24952568
E-mail: info@bionime.com
www.bionime.com

BIONIME USA CORPORATION

10865 Rancho Bernardo Road
Ste 100, San Diego, CA 92127, USA
Tel: +1 858 481- 8485
Fax: +1 858 481- 8472
E-mail: info.us@bionime.com

BIONIME (Shanghai) Co., Ltd.

Room 1110 Huawen International
Mansion No 999 Zhongshan West
Road Shanghai China
Tel: 86 21 6209 9336
Fax: 86 21 6208 3537
E-mail: info.cn@bionime.com

BIONIME GmbH

Tramstrasse 16
9442 Berneck, Switzerland
Tel: +41(0)71 722 98 40
Fax: +41(0)71 722 98 43
E-mail: info@bionime.ch
www.bionime.ch

BIONIME Australia Pty Limited

Level7, 60 York Street, Sydney NSW 2000
Postal address: PO BOX K602,
Haymarket NSW 1240
Tel: 02 9262 6900
Fax: 02 9262 6922
E-mail: info.au@bionime.com

BIONIME (ShenZhen) Co., Ltd.

No.708 Beauty AAA Building West Side
of Renmin Road Long Hua Street Baoan
District Shenzhen China
Tel: 86 755 82146630~1
Fax: 86 755 82146632
E-mail: info.cn@bionime.com

BIONIME
Blood Glucose Monitoring System

GM550

Auto Coding (Patent Pending)

Backlight

Noble Metal Electrode Strip



Rightest™ GM550

swissdesign 



About Bionime:

Established in April 2003, BIONIME is specialized in the business of biotechnology and medical testing. It integrates the top-tier expertise in medical science, chemistry, electronics and mechanism to commercialize the patented technology of Blood Glucose Monitoring System. BIONIME operates 3 subsidiaries in Europe, America and Asia. The products have been highly recognized by outstanding clinical tests and well-known journals. With the support from our customers, we are confident to establish the global competitive leadership.

Core Competency:

● Patented Strip Structure

- Solid plastic strip by injection
- Precise electrode assembly by insertion
- Leadless design securing signal measurement




Hydrophilic Cover film



Chemical Reagent
(GOD+Potassium ferricyanide+Polymer)



Au-BPE

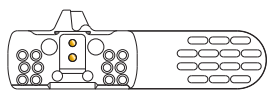


Injection Molding base (PBT)



Fabricated electrode Strip

Providing excellent base to reach high precision measurement.

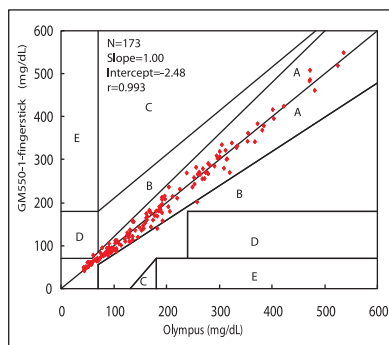

Source: Clinical Trial Report, Chung Shan Medical University Hospital, 2009

GM550 Meter:

The GM550 is a high-end model with many advanced functions. The feature of auto coding makes the meter easy to use. Backlight design enables people to read the test result even in the dark. With its accuracy & precision, the GM550 has attracted a lot of attentions from the professional users.



- Auto coding (Patent Pending)
- Backlight
- Nobel Metal Electrode Strip
- High Accuracy & Precision
- Large display – easy to read
- Tiny blood sample size
- Accurate measurement in few seconds

GM550 Specifications:

- Blood sample volume: 1.0μL
- Detection time: 5seconds
- Memory: 500 sets
- Test Range: 10-600 mg/dl (0.6-33.3mmol/l)
- Size: 90.6 x 46 x 16.5 mm
- Weight: 53 g (with batteries)

Easy to use



Friendly-large strip hand bar



Less pain-1.0 μL blood sample only, automatically drawn in



Rapidly-results less than 8 seconds



Safety-Keep from touching blood sample

GD500 Lancing Device:

The GD500 is an innovative device with a stylish look. With state of the art technology, the lancet not only helps reduce pain but also prevents contamination. The GD 500 is a sleek and practical partner for diabetic patients.

GD500 Acquires Taiwan "Certificate of Good Design Product"

- Less Pain
- Design Innovation
- Convenient
- Safety



GM550 Package Content:

- GM550 meter
- 10pcs of test Strips
- Control Solution
- GD500 Lancing Device
- Getting Starter Guide
- Test Strip Insert
- User Manual
- 10pcs of Lancet2
- 2 CR2032 Batteries
- Carrying Case
- Log Book
- Emergency Card

Easy to Identify

